


BETTER WATER LLC

1232 Reverse Osmosis Unit

Medical Device

Now With
Stainless
Steel
Frame!

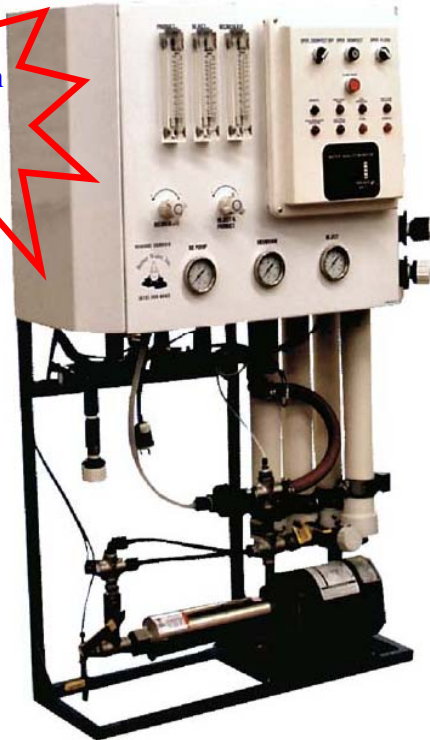

1232 RO UNIT

2800 Gallons per Day

The Better Water LLC perfect solution for acute units with 4-8 stations. This unit is built to the same high quality standards of our larger RO units and includes the same features. The 1232 series is economical, simple to operate, and quiet. The versatility of the unit allows for either tank or direct feed.

The 1232 RO carries a 510k certificate from the FDA and is designed to meet or exceed AAMI RD: 62 on RO quality water for dialysis. Each unit is individually tested prior to shipping.

Standard Features

- 5 micron filter (Pre-Filter)
- Water Quality Monitor
- Alarms
- Flow Meters
- Disinfect tank
- Product Water Divert-to-Drain
- Digital Displays
- Rugged Stainless Steel Frame Construction
- Baked on Powder Coat Finish
- Multi-staged air cooled stainless steel pump
- All RO units carry FDA510(k) #K920186C

Solutions. Service. Support.

Better Water LLC 698 Swan Dr Smyrna, Tn 37167 (615) 355-6063 and at www.betterwater.com

Equipment changes are at the sole discretion of BW,LLC. Brochure photographs & specs may not reflect recent changes



1232 R.O. Unit Medical Device

Models And Specifications

Model	EQRO-1232-2800
Capacity GPM/GPD	2.0/2880
# of Membranes	4
Size of Membranes	2.5"x 40" TFC
Control Voltage	24vac
Electrical Requirement	120vac/20 amp circuit
Dimensions	12"x32" x 56"
Incoming Water	4 gpm @ 30 psi
Feed Water Requirements	Temp:50f-90F,total chlorine <0.1ppm,SDI <3, TOC <2 ppm ,Operating pH range: 6-10,Bacterial Count <100 cfu
Gauges	RO Pump, Membranes, Reject
Flow meters	Product, Reject, Recirculation
Water Quality Monitor	Digitally Displays: % Rejection, Feed TDS, Product TDS, Alarm Set-Point
Standard Alarms	Poor Water Quality, Low Feed Pressure, High Pump Pressure, High Membrane Pressure, High Product Pressure and High Feed Water Pressure

BW LLC carries over 2000 parts, accessories and consumables in stock, and ready to ship including:

- Carbon Filters (provides 10 Minute Empty Bed Contact Time)
- Duplex & Simplex Softener
- Storage Tanks
- Brass pressure regulator
- Blend valve
- City boost pump
- Thermometer
- Floor valve boxes
- Remote alarms
- & More!



Central Water Systems	Central Delivery Systems	Reverse Osmosis Units	Portable RO Units
Pre-Treatment Racks	Post Treatment Racks	D.I. Exchange Accessories	Consumables
Hemodialysis	Medical	Pharmaceutical	Laboratories
		Commercial	Industrial